RACTIVEONE

RA-01 · 2026 EDITION

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
THE DEVICE

How it works

The Revive Active is an air compression leg massager — designed to relax and regenerate your leg muscles.

Four air chambers work in sequence. Thighs, calves, and ankles are massaged step by step — the technique of a manual massage, rebuilt in air pressure.

The result: improved blood circulation, relieved tension, real recovery for tired legs.

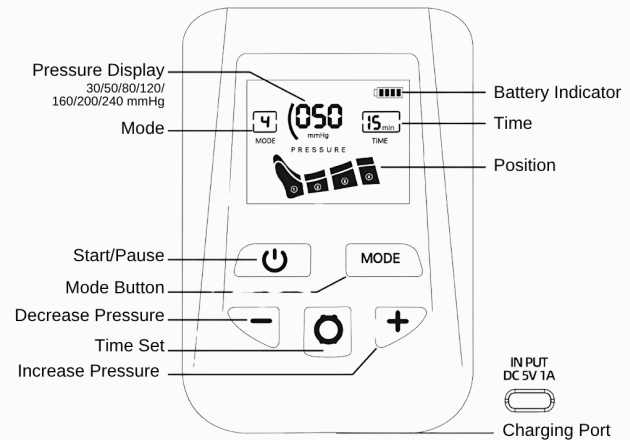START

How to start your recovery


- 1 Turn on**
Press the Start/Pause button to activate the controller.
- 2 Set pressure**
Use +/- to set the intensity. Start low, then work your way up.
- 3 Set duration**
Choose 15, 30, or 45 minutes via the Time button.
- 4 Find your mode**
Cycle through the programs until you find the right one.
- 5 Lean back**
The device takes it from here.


"The work is quiet. Most of it looks like this."

CONNECTION

How to connect the controller

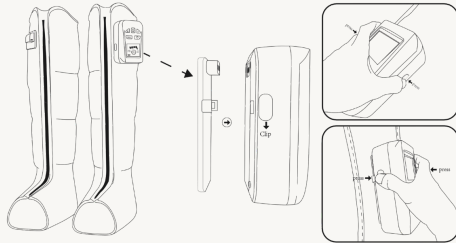Step 1 — Fix the controller. Place the controller on the bracket at the boot until it clicks in.

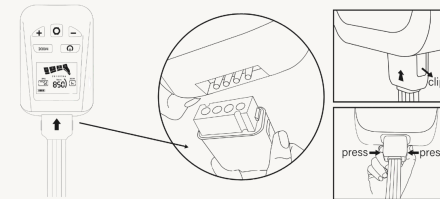Step 2 — Plug in the hose. Press the clips on both sides and push the hose in until it locks.

NOTICE AFTER USE

When storing, carefully release the hose connection — **always press the clips on both sides**. This prevents the connector from breaking.

MESSAGE MODES

Three programs, one goal

Each mode massages the air chambers in a different sequence. You choose what your body needs today.

MODE 01**Sequence**

Automatic, circular massage from bottom to top. The standard — ideal for daily recovery after training.

MODE 02**Circulation**

Individually selectable massage area. Direct the intensity to a spot that needs more work today.

MODE 03**Full Body**

The entire leg is massaged step by step from bottom to top — for a full activation of blood flow.

PRESSURE LEVELS**30 · 50 · 80 · 120 · 160 · 200 · 240 mmHg**

BATTERY & CARE

Indicators and daily care

BATTERY LEVEL

● Green	100 – 65 %
● Blue	65 – 25 %
● Red	25 – 5 %

Below 5 % the device shuts off automatically. While charging, the display shows "CH" — the device cannot be used during charging.

CLEANING & STORAGE

- ✓ Wipe the controller with a soft cloth.
- ✓ Clean fabric with alcohol or a damp cloth.
- ✗ No detergents, no immersion cleaning.
- ✓ Store dry — out of reach of children.
- ✓ Recharge every two months when not in use.

TROUBLESHOOTING

When something isn't working

DISPLAY / PROBLEM	POSSIBLE CAUSE	SOLUTION
Device won't start	Battery is empty.	Charge the device before turning it on.
Device won't start	Timer active — the device has shut off automatically.	Restart the device.
Er1 on display	Overvoltage alarm.	After automatic shutdown and voltage reduction, restart.
Low pressure value	Air leak in the product.	Contact customer service within the warranty period.

NEED HELP?

Write to us at support@ractiveone.com — we respond personally on weekdays, usually within 24 hours.

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